


Cataract-Bi-manual Micro phaco Set

DETAILS OF Akahoshi 3-prong Speculum

Catalog No	AE-1023
Name	Akahoshi 3-prong Speculum
Category	Specula
Description	3 prong, Nasal speculum for easy temporal access, adjustable self locking mechanism. Solid prong: 38mm, Open prong: 42mm.

Image of Akahoshi 3-prong Speculum

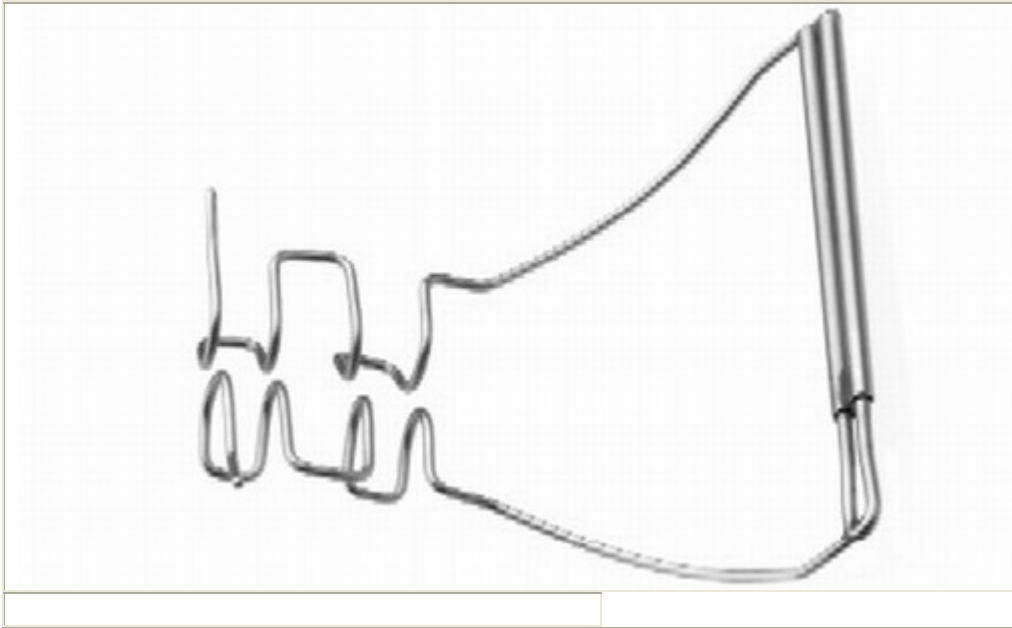

DETAILS OF Ikeda II Micro Capsulorrhexis Forceps

Catalog No	AE-4385
Name	Ikeda II Micro Capsulorrhexis Forceps
Category	Capsulorrhexis Forceps
Description	23 gauge (.060mm), Overall Length:139mm. Sharp Tip can be used as Cystotome. 90 degree curved forceps allow for easy maneuverability through the paracentesis incision. Specially designed platform in the jaws allow for a secure hold of the anterior capsule in the presence of viscoelastic. Round handle provides better control.

Image of Ikeda II Micro Capsulorrhexis Forceps

