


## Diagnostic Lenses



### Single Mirror Gonio Lens

Small size gonio lens with one 62° mirror. Compact knurled ring simplifies 360° viewing of the anterior chamber angle.

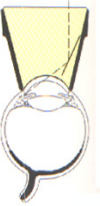

Product Code	Contact Diameter	Lens Height	Lens Magnification
OSMG	15mm	19.3mm	.93x



### Two Mirror Gonio Lens

Two opposing 62° mirrors provide a complete view of the anterior chamber angle with only a 180° lens rotation. Methylcellulose and no methylcellulose designs available.




Product Code	Contact Diameter	Lens Height	Lens Magnification
O2M	15mm	19.3mm	.93x
O2MG*	15mm	19.3mm	.93x



### Four Mirror Mini Gonio Lens

Four 62° mirrors allow complete observation of the angle with little lens rotation. Small diameter flange is convenient for eyes with small palpebral fissures. Anterior holding ring available in small and large sizes.

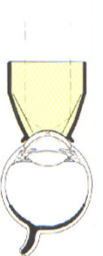


Product Code	Contact Diameter	Lens Height	Ring Diameter
O4GF	15mm	22.1mm	23mm
O4GF-LR	15mm	26.6mm	32mm



### Thorpe Four Mirror Gonio Lens

Four 62° mirrors give a 360° view of the anterior chamber angle with only a slight lens rotation. Posterior pole can be viewed through center of lens.

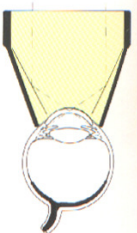

Product Code	Contact Diameter	Lens Height
OT4MG	18mm	32.2mm

\*no methylcellulose required